
- Ⓓ **Montage- und Gebrauchsanleitung**
- ⒼⒷ **Installation and Operating Instructions**
- Ⓕ **Notice de montage et mode d'emploi**
- ⒶⒻ **Installatie-en gebruikshandleiding**
- Ⓘ **Istruzioni di montaggio e per l'uso**
- Ⓔ **Instrucciones de instalación y funcionamiento**
- ⒶⒻ **Instrukcja montażu i obsługi**

**D Inhalt**

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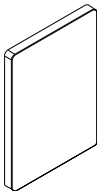### **Purpose according to Directive 93/42/EEC**

The HEWI support rail is designed for use in bathrooms and lavatories. It is used as an aid to sitting down and standing up in lavatories and also offers people with physical disabilities additional support and something to hold on to.

### **Warning and hazard notices**

- The HEWI support rail is designed for a maximum static load of 100 kg vertically. After fitting the support rail, make sure that the lock is properly engaged.
- Improper mounting or use can lead to accidents involving serious injuries.
- A certain amount of vertical and horizontal bending is normal when the support rail is under load. In the event of damage (e.g. permanent bending), the support rail may only continue to be used following inspection by a specialist company responsible for maintenance.
- Do not misuse the support rail, e.g. as gymnastics apparatus or climbing aid.
- When folding down the support rail, ensure that nobody is located beneath the rail.

### **Care and cleaning instructions**

HEWI products are easy-care. Generally speaking, an occasional wipe with a damp cloth is enough. When using detergents, please observe the following:

- Only use detergents with a pH of between 6 and 8.
- Detergents containing active substances such as phosphates, soap and tensides may be used.
- Follow the detergent manufacturer's instructions (concentration, contact time).
- Detergents containing acids, alkalis, bleach or scouring agents may not be used.
- Do not use scouring cleaning equipment, such as brushes

### **Exclusion of liability**

HEWI accepts no liability whatsoever for damage to persons or property that occurs or has occurred as a result of improper mounting or maintenance, particularly if the safety instructions and other points to note listed are not or have not been observed in full. HEWI accepts no liability whatsoever for damage occurring from incorrect cleaning.

**For further information, please contact HEWI Technical Support,  
tel. +44 1634 377688.**

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HEWI Heinrich Wilke GmbH